# **ION Health Suite**

## **SOFTWARE MANUAL**

∴ ENGLISH ( 3 – 14 ) ∴

## **MANUAL DE SOFTWARE**

∴ ESPAÑOL ( 15 – 26 ) ∴

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∴ FRANÇAIS ( 27 – 38 ) ∴

## **MANUALE PER L'USO DEL SOFTWARE**

∴ ITALIANO ( 39 – 50 ) ∴

## **SOFTWAREHANDBUCH**

∴ DEUTSCH ( 51 – 62 ) ∴

## **SOFTWAREHANDLEIDING**

∴ NEDERLANDS ( 63 – 74 ) ∴

MANUAL VERSION 1.0  
2011.08.31

## INSTALLING THE SOFTWARE


### MINIMUM SYSTEM REQUIREMENTS

Before starting the installation, make sure your computer meets (or exceeds) the following minimum system requirements to ensure the software runs properly:

- PC with Microsoft® Windows 98, 2000, XP, Vista (32-/64-bit), or 7 (32-/64-bit)
- 400 MHz processor (or higher)
- 64 MB RAM (or more)
- 75 MB free hard drive space
- 1024 x 768 graphics card (or higher resolution)
- Available USB port (USB 1.1 or 2.0)
- **Optional:** Windows-compatible printer (for printing copies of records)

### SOFTWARE INSTALLATION

1. Insert the included CD into your computer's CD drive. The installation window should appear automatically. Click **Install ION Health Suite Software**.
  - If it does not appear, double-click **My Computer** and click on your computer's CD drive.
  - If the installation window does not appear after that, right-click the CD drive, select **Explore**, and double-click the **Setup** file.
2. Follow the on-screen instructions to install the software.
3. When you are finished, click **Close** to close the installer.
4. Open ION Health Suite by double-clicking the shortcut icon on your Desktop or by going to **Start ► All Programs ► ION HEALTH SUITE ► ION HEALTH SUITE**.
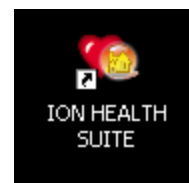
## **IMPORTANT NOTES**

### **About ION Health Suite**

ION Health Suite is a personal computer (PC) software application that lets you download readings from your ION USB Wrist Blood Pressure Monitor or USB Insta-Scan Thermometer, store them securely on your computer, and analyze them with graphs and other software tools. This information can provide you with a better understanding of your blood pressure or temperature to help you and your health care professional manage your health more effectively.

### **ION Health Suite allows you to:**

- Transmit blood pressure or temperature results from a device to your PC
- Track your readings digitally in the software instead of in a log book
- Present the trends and statistics of your readings
- Store your readings securely on your PC
- Print and share your data with your doctors

### **Warning:**

- The readings taken by the ION USB Wrist Blood Pressure Monitor or USB Insta-Scan Thermometer are for reference only. Please consult your physician to discuss your results before adjusting your diet, exercise program, etc.
- To avoid the possibility of electrical shock, never perform a blood pressure or temperature test while the device is connected to the computer.
- When using ION Health Suite, make sure you have selected the correct account. ION Health Suite may not properly download one user's results to another user's account.

## GET STARTED

### CONNECT YOUR DEVICE

Before downloading any readings to your computer or setting up the blood pressure monitor or thermometer, follow these steps:

1. With the blood pressure monitor or thermometer powered off, connect it to a free USB port on your computer using the included mini-USB cable.
2. "PCL" ("PC link") will appear in the device's display, indicating it is connected properly.
3. Open ION Health Suite by double-clicking the shortcut icon on your Desktop or by going to **Start ► All Programs ► ION HEALTH SUITE ► ION HEALTH SUITE**.

### CREATE A USER PROFILE


1. Enter your name (up to 50 characters), birth date, and sex. (After clicking the Birthday drop-down menu, you can jump to a different month or year by clicking on them at the top of the calendar.) Click the ➡ button to continue.
2. In the window that appears, click the drop-down menus to select the desired setting:
  - **Language:** English, Spanish, French, Italian, German, or Dutch
  - **Height Unit:** feet/inch (U.S.) or centimeter (metric)
  - **Weight Unit:** pound (U.S.) or kilogram (metric)
  - **BP (Blood Pressure) Unit:** mmHg (recommended) or kPa
  - **Thermometer Unit:** F° or C°
  - **Device:** BP Monitor (for USB Wrist Blood Pressure Monitor) or Thermometer (for USB Insta-Scan Thermometer)
3. Click the ➡ button to save the data and continue to the main window in the software.

#### Notes:

- To cancel the setup at any time, click the "X" (Cancel) button in the lower right corner of the window.
- You can change your name later. See SET YOUR PERSONAL SETTING in the ADJUST YOUR SETTINGS section for more information.
- This setup window only appears automatically the first time you use the software. Once a user has been created, the software will open to the main window in the future. To create more users, click "**Add new users**" on the main window.
- You can enter a maximum of five users. To delete users, click "**Delete users**" on the main window, check the box next to the user(s) you want to delete, and click the trash can button to confirm deletion, or click the "X" button in the lower right corner of the window to cancel and return to the main window.

## SOFTWARE OVERVIEW

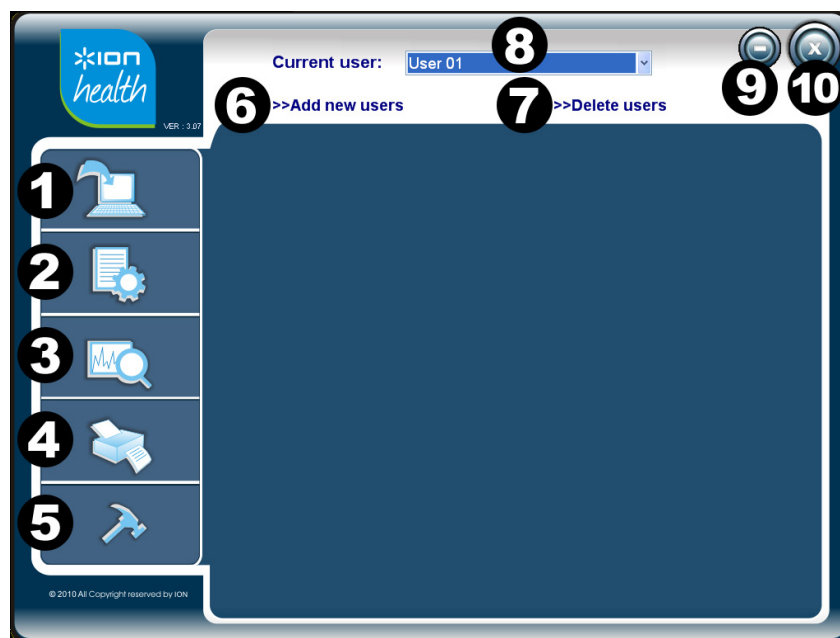1. **Download Menu** – Place your mouse over this menu and click **Download** to download data from your device.
2. **Data Menu** – Place your mouse over this menu and click **BP Data** (blood pressure) or **Thermometer Data** (temperature) to review your data in a table. See the VIEW YOUR DATA IN A TABLE section for more information.
3. **Graph Menu** – Place your mouse over this menu and click **BP Graph** (blood pressure) or **Thermometer Graph** (temperature) to review your data in a graph. See the VIEW YOUR DATA IN A GRAPH section for more information.
4. **Print Menu** – Place your mouse over this menu and click **Print Data (BP)** (blood pressure) or **Print Data (Thermo.)** (temperature) to print a report of your data. See the PRINT YOUR DATA section for more information.
5. **Setting Menu** – Place your mouse over this menu and click one of the following options (see the ADJUST YOUR SETTINGS section for more information):
  - **Personal Setting:** Edit your user information.
  - **Device Setting:** Adjust the settings for your device (blood pressure monitor or thermometer).
  - **Software Setting** Adjust the settings for ION Health Suite.
6. **Add new users** – Click this to create more user profiles. You will see the same window as you did during the software's initial setup. See the CREATE A USER PROFILE for more information.
7. **Delete users** – Click this to delete a user.
8. **Current user** – Click this drop-down menu to switch to another user profile. (The software can store a maximum of 5 users.)
9. **Minimize button** – Click this to "minimize" the window. Click the "ION HEALTH SUITE" bar in your computer's Taskbar to restore the window.
10. **Close button** – Click this to close the software.

## ADJUST YOUR SETTINGS

Follow the steps in this section to customize the settings of:

1. Your user information (see **SET YOUR PERSONAL SETTING**)
2. Your ION USB Wrist Blood Pressure Monitor or USB Insta-Scan Thermometer (see **SET YOUR DEVICE SETTING**)
3. The software itself (see **SET YOUR SOFTWARE SETTING**)

### SET YOUR PERSONAL SETTING

This window contains two tabs: **Personal Information** and **Measurement Interval**. These parameters will be reflected when you review your readings.

1. Place your mouse over the **Settings Menu** (hammer icon) and click **Personal Setting**.
2. The first tab is the **Personal Information** tab.

- Enter your name, birthday, gender, height, and weight.
- **If you are using the ION USB Wrist Blood Pressure Monitor:** The fields for **Sys Pressure (Lower/Upper Limit)** and **Dia Pressure (Lower/Upper Limit)** are the recommended minimum / maximum values for your systolic and diastolic pressures (respectively). The values you first see in these fields have been set according to recommendations from the American Diabetes Association. However, we highly recommend consulting with your health care professional to determine if these values are ideal for you.

- **If you are using the ION USB Insta-Scan Thermometer:** The fields for **Thermometer Level (Lower/Upper Limit)** are the typical minimum/maximum values for the temperature of the human body. We highly recommend consulting with your health care professional to determine if these values are ideal for you.

3. Click the **Save** button (disk icon) in the lower right corner to save your information.

4. Click the **Measurement Interval** tab.

5. Each row represents an interval of the day during which you plan to take a blood pressure reading or temperature reading. Click the "AM/PM," the hour, and the minute in the white fields, then type in the number or use the up/down arrows to set the end time of that interval. The gray fields will adjust automatically.

- **Example:** If you plan to take a reading before breakfast, which you normally do at 8 a.m., set the "Before breakfast" interval to end at 8 a.m.
- **Note:** These intervals cannot overlap each other or be set out of sequence. If you try to set the end time of an interval to overlap another, you will see an error message. Simply click "OK" and set the interval so it does not overlap another.
- **Note:** You can click the Default button (paper/star icon) to return all the intervals to their default times.

6. Click the **Save** button (disk icon) in the lower right corner to save your information.

## SET YOUR DEVICE SETTING

This window allows you to do two things:

1. Delete all readings from the device's memory
2. Set the device's time.

These are independent functions. If the time is reset, the readings currently stored in the device will not be changed; they will retain their old time settings. Similarly, erasing the device's memory will not affect its current date or time.

1. With the device turned off, use the included mini-USB cable to connect it to a free USB port on your computer.
2. Place your mouse over the **Settings Menu** (hammer icon) and click **Device Setting**.
3. The first tab is the **Delete the device's memory** tab. (To skip this and set the device's time instead, skip to Step 7 below.)
4. Click the **Device** drop-down menu and select your device – **BP Monitor** or **Thermometer**.
5. Click the ➡ button. A confirmation window will appear if the connection is successful. Click **OK** to continue.
6. Click the ➡ button. A second confirmation window will appear, asking if you are sure you want to delete the device's readings. Click **Yes** to continue (or **No** to cancel).
7. To set the device's time, click the **Set the device's time** tab.
8. Click the **Device** drop-down menu and select your device – **BP Monitor** or **Thermometer**.
9. Click the ➡ button. A confirmation window will appear if the connection is successful. Click "OK" to continue.
10. Click one of the two boxes shown in the window:
  - **Time:** The device will have the same time setting as your computer.
  - **Time setting:** Set the device's time manually by clicking the **AM/PM**, the hour, and the minute in the white fields, then type in the number or use the up/down arrows.
11. Click the ➡ button.

**Note:** ION Health Suite does not account for Daylight Saving Time, so you will need to adjust your device's time accordingly.

## SET YOUR SOFTWARE SETTING

This window allows you to set what language and measurement units are shown in the software.



1. Place your mouse over the **Settings Menu** (hammer icon) and click **Software Setting**.
2. Click the drop-down menu next to each item to set it:
  - **Language:** English, Spanish, French, Italian, German, or Dutch
  - **Height Unit:** feet/inch (U.S.) or centimeter (metric)
  - **Weight Unit:** pound (U.S.) or kilogram (metric)
  - **BP Unit:** mmHg (recommended) or Kpa
  - **Thermometer Unit:** F° or C°
  - **Device:** BP Monitor (for USB Wrist Blood Pressure Monitor) or Thermometer (for USB Insta-Scan Thermometer)
  - **User No.:** n/a
3. Click the **Save** button (disk icon) in the lower right corner to save your information.

## DOWNLOAD YOUR DATA

To download readings from your device to your computer, follow these steps:

1. Place your mouse over the **Download Menu** (computer icon) and click **Download**.
2. With your device powered off, use the included mini-USB cable to connect it to a free USB port on your computer.
3. Click the **Device** drop-down menu and select your device – **BP Monitor** or **Thermometer**.
4. Click the ➡ button.
5. Click the drop-down menu and select how many readings you want to download from your device to your computer.
6. Click the **Download** arrow button (in the lower right corner) to start downloading the data, or click the ⬅ button to return to the previous window.
7. After all the data is downloaded, it will be shown on the window:

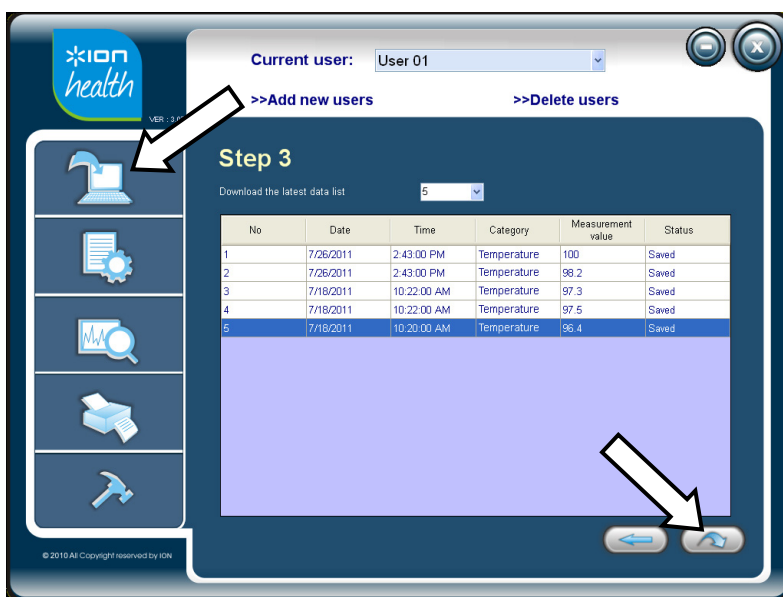
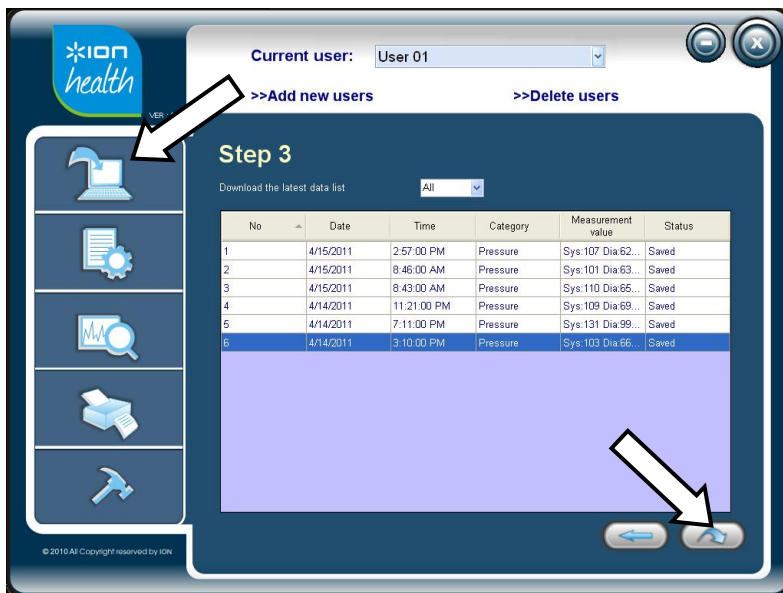
- **No (Number):** The number of the reading in the sequence.
- **Date:** The date the reading was taken (according to the device's date setting).
- **Time:** The time the reading was taken (according to the device's time setting).
- **Category:** If the reading was taken from a blood pressure monitor, the **Category** will be **Pressure**. If the reading was taken from a thermometer, the **Category** will be **Temperature**.
- **Measurement Value:** If the reading was taken from a blood pressure monitor, the reading will include your blood pressure values (systolic pressure, diastolic pressure, and pulse rate). If the reading was taken from a thermometer, the value will be your temperature.
- **Status:** If this is the first time the reading has been downloaded to your computer, it will say **Saved**. If it has already been stored on your computer, it will say **Exist**.

8. When you are finished, you can disconnect your device from your computer. The readings you downloaded can now be viewed in the software.

## VIEW YOUR DATA IN A TABLE

To view your readings in a table, follow these steps:

1. Place your mouse over the **Data Menu** (paper/gear icon) and select which data you want to view in a table – **BP Data** (blood pressure) or **Thermometer** (temperature).
2. Check the box next to **Select the date interval** or **Select a specific date range**. Click the drop-down menu(s) next to it to select the range of dates whose readings you want to view. (When setting a specific date, you can jump to a different month or year by clicking on them at the top of the calendar.)
3. Check the box next to **Select the time interval** or **Select a specific time range**. Click the drop-down menu(s) next to it to select the range of times whose readings you want to view.
4. Click **OK** to view your readings, which will appear in the table below.
  - **Date / Time:** The date and time the reading was downloaded to your computer from your device.
  - **Systolic Pressure / Diastolic Pressure / Pulse or Thermometer:** Your blood pressure or temperature reading.
  - **Note:** You can double-click this column and a new window will appear, allowing you to type in any additional information you would like to store with this reading.
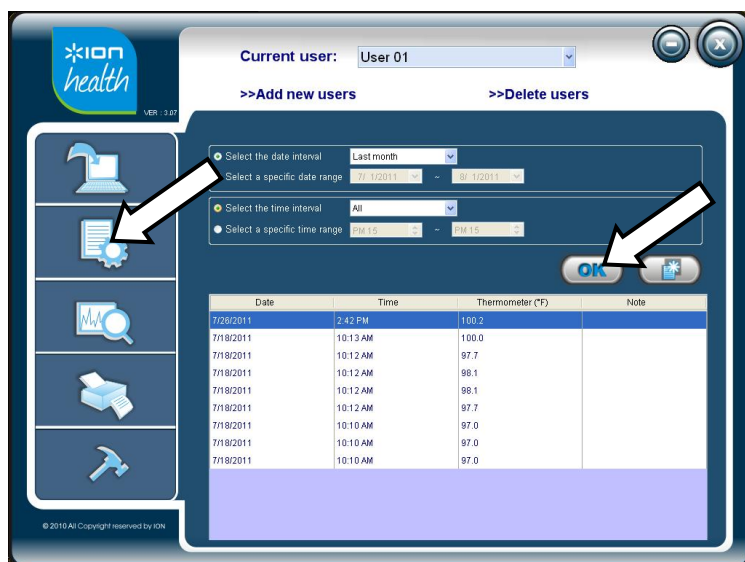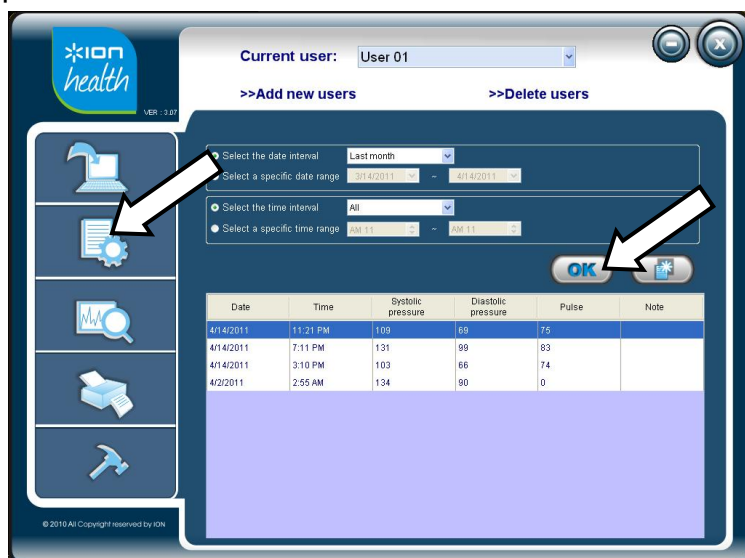
## TO EXPORT DATA (OPTIONAL)



1. Click the "paper/star" icon (next to **OK**).
2. In the window that appears, click the large white field under **Export File Name**.
3. You will be asked to select a location to save the file and the file name. (We recommend saving it to someplace obvious like your Desktop and naming it with the date so you can keep track of it easily.) When you have done so, click **Save**.
4. To access the file with your data, go to that location on your computer. You can use Microsoft Excel to open the .csv file.

## VIEW YOUR DATA IN A GRAPH

To view your readings in a graph, follow these steps:


1. Place your mouse over the **Graph Menu** (graph/magnifying glass icon) and select which data you want to view in a graph – **BP Data** (blood pressure) or **Thermometer** (temperature).
2. Check the box next to **Select the date interval** or **Select a specific date range**. Click the drop-down menu(s) next to it to select the range of dates whose readings you want to graph. (When setting a specific date, you can jump to a different month or year by clicking on them at the top of the calendar.)
3. Check the box next to **Select the time interval** or **Select a specific time range**. Click the drop-down menu(s) next to it to select the range of times whose readings you want to graph.
4. Click the **Grouping Unit** drop-down menu and select how you want the readings to be grouped on the graph (by **Date**, by **Time**, or by **Day** [of the week]).
5. Click **OK** to view your readings, which will appear in the table below.
  - **X-axis (Date, Time, or Day):** This is the interval you selected as the **Grouping Unit**.
  - **Y-axis (Pressure or Thermometer):** Your blood pressure or temperature reading (respectively).
  - **Data Points:** You can place your mouse over any of these points on the graph to view its specific information.
    - **Blood Pressure:**
      - **Blue dots** represent individual systolic pressure readings.
      - **Red diamonds** represent the average of all systolic pressure readings in that column.
      - **Green triangles** represent individual diastolic pressure readings.
      - **Purple squares** represent the average of all diastolic pressure readings in that column.
    - **Temperature:**
      - **Blue dots** represent individual temperature readings.
      - **Red squares** represent the average of all temperature readings in that column.
  - **Target Range Area:** The shaded band across the graph indicates your lower and upper limits of your blood pressure or temperature (which you can set in the Personal Settings). (For temperature, the shaded band will be gray. For blood pressure, the systolic pressure range appears in orange while the diastolic pressure range appears in gray.)
  - **Out-of-Target Area:** This yellow-shaded area represents the range beyond the lower and upper limits of your blood pressure or temperature (which you can set in the Personal Settings).

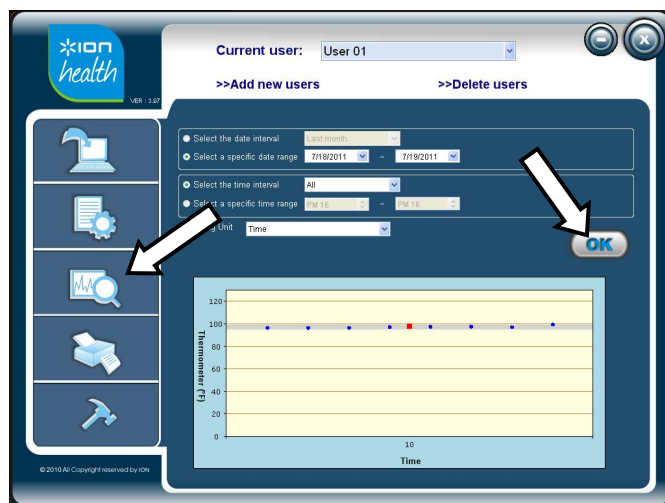
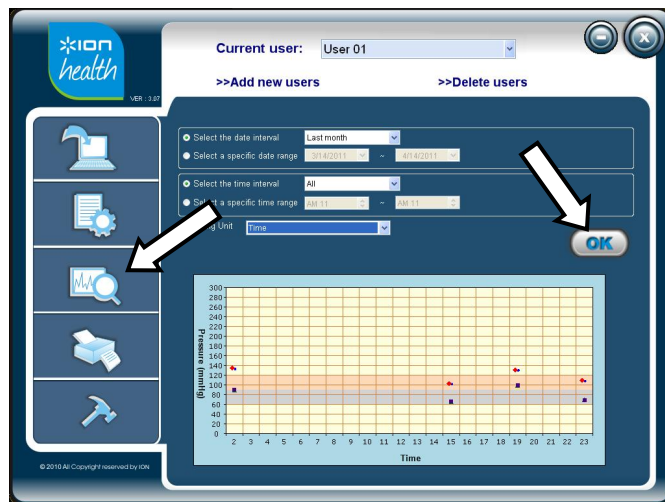
## PRINT YOUR DATA

To print your readings, follow these steps:

1. Place your mouse over the **Print Menu** (printer icon). To print blood pressure data, click **Print Data (BP)**. To print temperature data, click **Print Data (Thermo.)**.
2. Check the box next to **Select the date interval** or **Select a specific date range**. Click the drop-down menu(s) next to it to select the range of dates whose readings you want to print. (When setting a specific date, you can jump to a different month or year by clicking on them at the top of the calendar.)
3. Check the box next to **Select the time interval** or **Select a specific time range**. Click the drop-down menu(s) next to it to select the range of times whose readings you want to print.
4. Skip the **Report Name** drop-down menu. (Your type of data – **Blood pressure report** or **Thermometer report** – will be the only option.)
5. Click the **Grouping Unit** drop-down menu and select how you want the readings to be grouped on the graph (by **Date**, by **Time**, or by **Day** [of the week]).
6. Click **OK**.
7. Your data will be compiled in a report with your user information (name, height, weight, etc.) and all your readings for the designated date/time range. The data will be shown in table and graph forms in the report (on multiple pages).

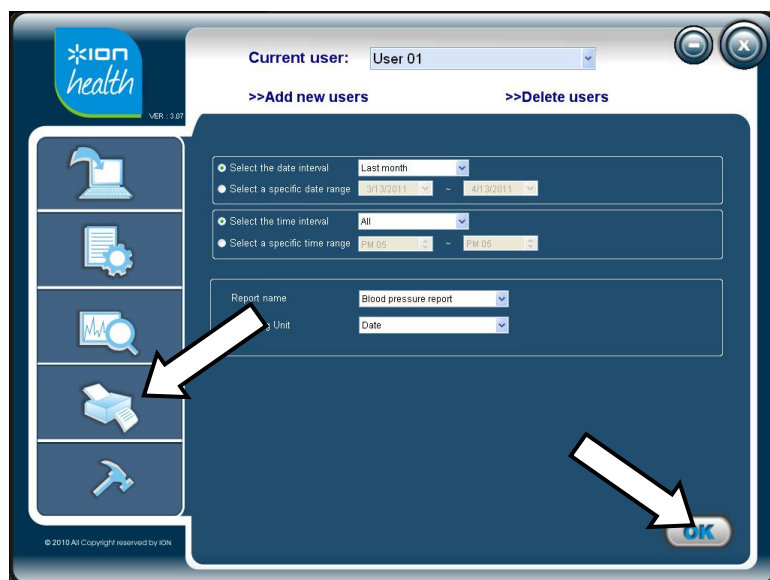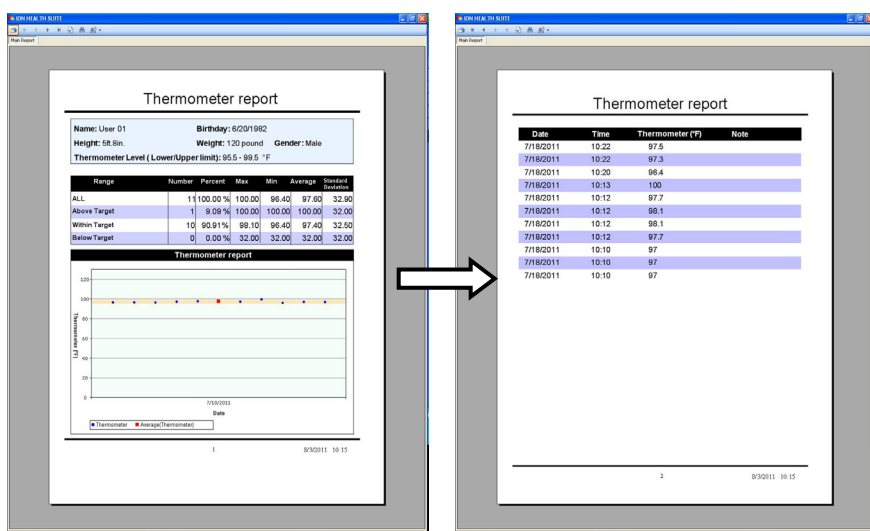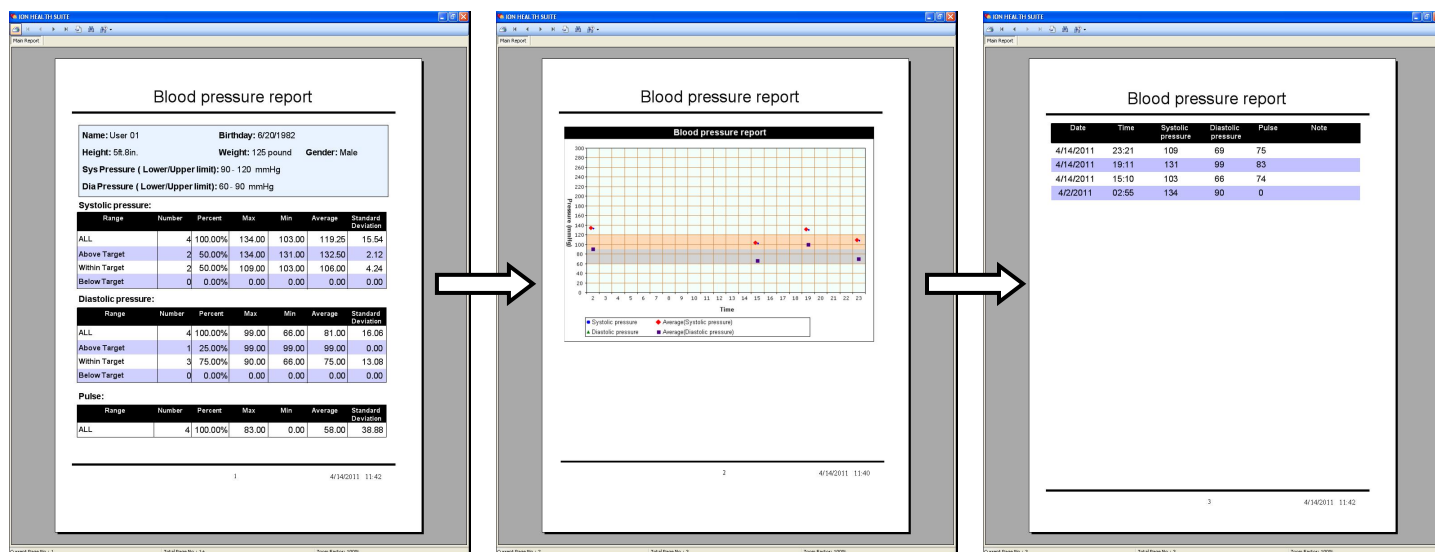1 3 2 2 3 4 5 6

Use the buttons in the toolbar at the top of the window to do the following:

1. **Print:** Click this to print the report. (Your computer must be connected to a printer that is powered on.)
2. **Go to Previous/Next Page:** Click these to move to the previous/next page.
3. **Go to First/Last Page:** Click these to skip to the first or last page.
4. **Go to Page:** Click this and enter a page number to skip to that page.
5. **Find Text:** Click this and enter a word to skip to the page(s) containing that word.
6. **Zoom:** Click this and select how much you want to zoom in/out of the page.

## Notes:

- Do not make any adjustments to your diet, exercise program, medication, etc. based on this data – consult your health care professional first.
- This report is only for printing and cannot be saved directly to your computer. To save your data, see the "To export your data" instructions in the VIEW YOUR DATA IN A TABLE section.

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